


Inclusive Education

AAC App Set up for Touch Chat

iPads should initially be locked into Landscape orientation

How:

- Draw down from top right corner to access the control panel, press on the lock once you have it in the landscape orientation.

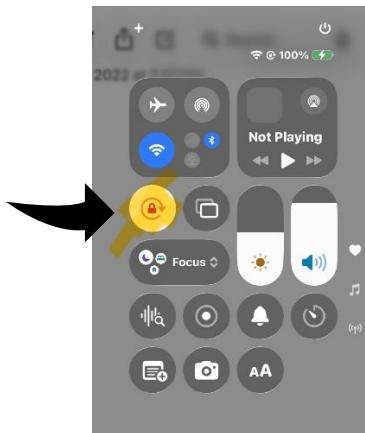

When starting with a brand-new Vocabulary:

- To choose a new Vocabulary, go to the Vocab section located at the top left corner of the page then decided which WORD POWER vocabulary file to use (typically this is WordPower60SS)
- Customize this vocabulary and rename it 0 – the student's first name and then last name Initial – Month2025 (e.g. 0 - KrisM – Dec2025)

Touch Chat Settings (accessed through the Menu in the app):

<p>Done Settings</p> <p>Profiles</p> <p>Save Profile ></p> <p>Load Profile ></p> <p>Navigation</p> <p>Page Animation <input checked="" type="checkbox"/></p> <p>Allow Gestures <input checked="" type="checkbox"/></p> <p>Gesture Pop Sound <input checked="" type="checkbox"/></p> <p>Button Click Sound <input checked="" type="checkbox"/></p> <p>Button Navigation Icons <input checked="" type="checkbox"/></p> <p>Page Size</p> <p>Status Bar <input checked="" type="checkbox"/></p> <p>Keyguard Inset Off ></p>	<p>Done Settings</p> <p>Editing</p> <p>Menu Passcode</p> <p>Use Menu Passcode <input type="checkbox"/></p> <p>Editing Note Editing Note: Menus may be completely disabled in the Apps section of the iOS Settings. (Allow Editing)</p> <p>Allow Vocabulary Deletion <input checked="" type="checkbox"/></p> <p>Speech</p> <p>English voice Melany ></p> <p>French (Canada) voice Amélie ></p> <p>Phone/FaceTime Voice Samantha ></p> <p>Novelty Voices <input checked="" type="checkbox"/></p> <p>Speech Display Bar</p>	<p>Options for Voices:</p> <ul style="list-style-type: none"> • Julian (boy) • Victoria (girl) • Melanie (woman) • Elijah (teen male) US
---	--	---

Done	Settings
Speech Display Bar	
Hide Speech Display	<input checked="" type="checkbox"/>
Allow Copy/Paste	<input checked="" type="checkbox"/>
Auto Capitalization	<input checked="" type="checkbox"/>
Show Icons	<input checked="" type="checkbox"/>
Hide Delete Button	<input checked="" type="checkbox"/>
Font	Arial >
Font Size	50 Point >
Font Color	Black >
Background Color	White >
Number of Lines	1 Line >
Visible Cursor	<input checked="" type="checkbox"/>

Done	Settings
Speech Display Bar	
Number of Lines	1 Line >
Visible Cursor	<input checked="" type="checkbox"/>
Allow Place Cursor	<input checked="" type="checkbox"/>
Tap Action	Speak >
Highlight Spoken Word	OFF >
Highlight Color Phone/FaceTime	
Expanded Speech Area	
Font	Arial >
Font Size	80 Point >
Font Color	Black >
Background Color	White >

Done	Settings
Speak	
Speak Characters	<input checked="" type="checkbox"/>
Speak Words	<input checked="" type="checkbox"/>
Speak Sentences	<input checked="" type="checkbox"/>
Auto Clear	Off >
Grammar	<input checked="" type="checkbox"/>
Prevent Interruption	<input checked="" type="checkbox"/>
Data Logging	
Enable Logging	<input checked="" type="checkbox"/>
Privacy Password	
Realize Language™ Setup	
Buttons	
Activate on Release	<input checked="" type="checkbox"/>

Done	Settings
Buttons	
Activate on Release	<input checked="" type="checkbox"/>
Dwell Time	None >
Release Time	None >
Force Font Size	<input checked="" type="checkbox"/>
Highlight Touch	<input checked="" type="checkbox"/>
Gradient Fill	<input checked="" type="checkbox"/>
Gradient Saturation	<input checked="" type="checkbox"/>
Hide All Icons	<input checked="" type="checkbox"/>
Show Hidden Buttons	<input checked="" type="checkbox"/>
Show grid	<input checked="" type="checkbox"/>
Dynamic Grammar Labels	<input checked="" type="checkbox"/>
Button Margin	5 Point >

Done	Settings
Prediction	
Prediction Type	Learning + Static >
Learned Words	
Whiteboard	
Disable drawing in split view	<input checked="" type="checkbox"/>
Scale drawing in split view	<input checked="" type="checkbox"/>
Reload last drawing	<input checked="" type="checkbox"/>
Autosave drawing	<input checked="" type="checkbox"/>
Head Tracking	
Enable Head Tracking	<input checked="" type="checkbox"/>
Social Networking	
Allow Social Networking	<input checked="" type="checkbox"/>
Amp Connect	
Select Amp	Not Connected >